


Axessio™

Peripheral microguidewire

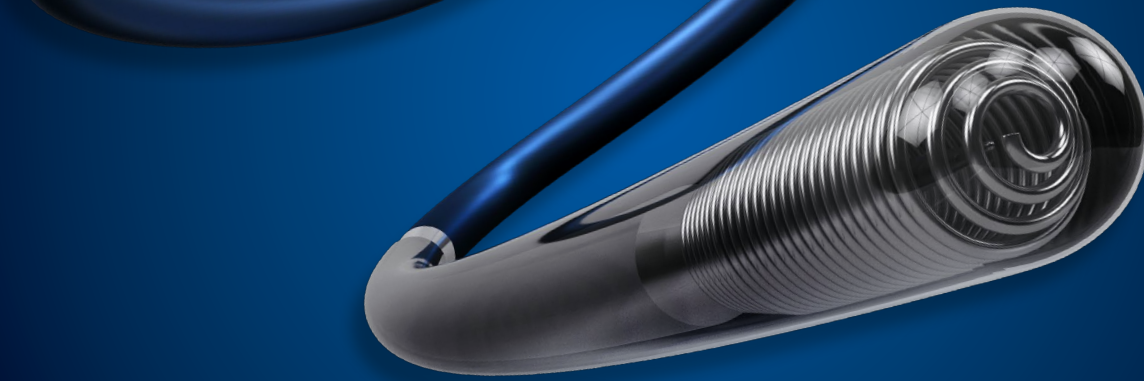

Go the distance. Access the result.

Axessio™ shapeable peripheral microguidewire enables navigation in challenging anatomies. Twister™ core technology is designed to provide predictable torque control for durable performance in complicated cases*.

CLINICAL APPLICATION VERSATILITY


- Shapeable tip to tailor to various anatomies
- Navigate challenging and highly tortuous anatomy maintaining precise 1:1 torque response

GREAT DURABILITY & RESISTANCE IN CHALLENGING CASES

- Resistance to torque deformation
- Stable torque control during long procedures

EXCELLENT VISIBILITY

- 39 cm highly radiopaque jacket
- 3 cm platinum coil tip

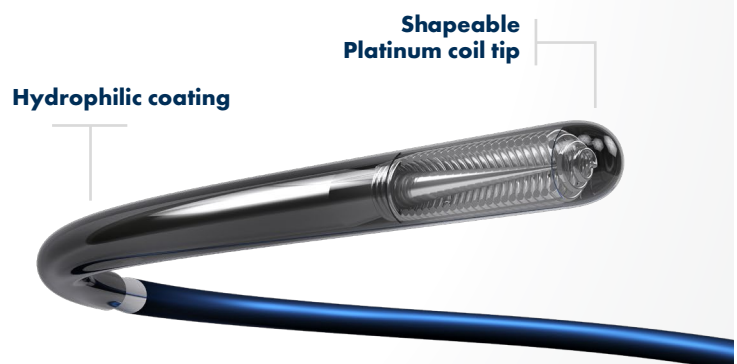


Equilibrium steerability & support

- Blend of flexibility and stiffness
- Core to tip design/1 piece construction
- Stainless steel core
- Twister™ core technology

* Integer internal data

Technical Specifications

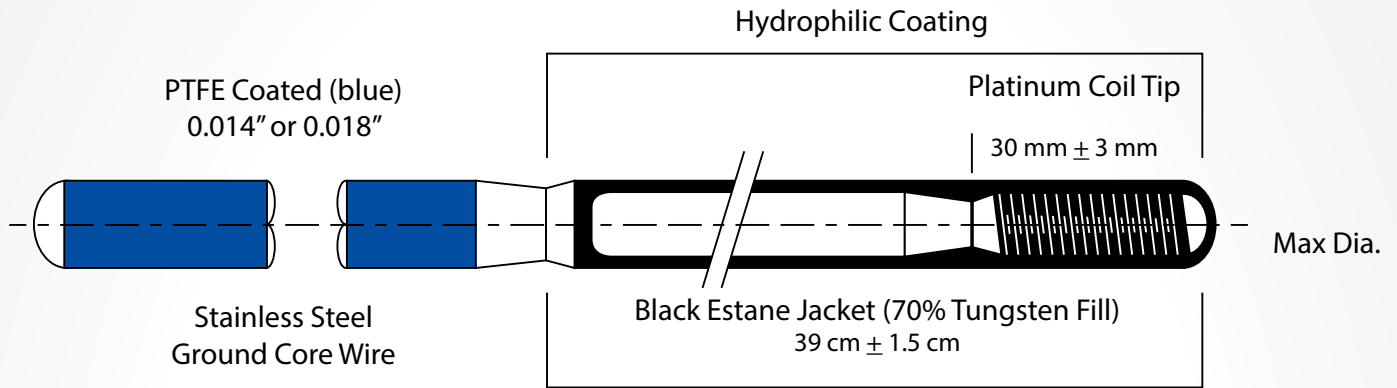

Distal O.D. (in)	Order code	Total length	Platinum Coil Tip Length (cm)	Radiopaque length (cm)	Tip Shape
0.014	240428	180cm	3cm ± 0.3cm	39 cm ± 1.5 cm	Straight
0.018	240429	180cm	3cm ± 0.3cm	39 cm ± 1.5 cm	Straight

Guerbet | 

interventional.guerbet.com